



# Danumed Overview Guide

*danumed button-sets*

*danumed extension set*

*danumed gastrotube*

*danumed stomaMeter*

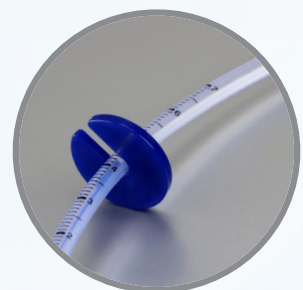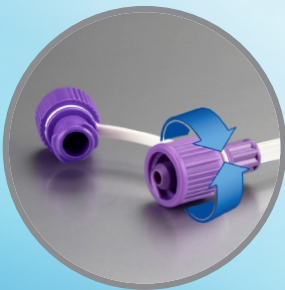

# danuButton®

## Compact Set / Premium Set

*The danuButton® is a replacement feeding tube for long-term enteral nutrition and use with completely healed stoma. The low-profile feeding tube is easy to conceal under clothing and especially for active, mobile patients. The extensive accessories of the danuButton® Sets enable you to an easy and comfortable replacement of the gastrostomy probe, additional components are no longer necessary.*

*The decision for the danuButton® is also the decision for more quality of life through an innovative, high-quality product.*



### COMFORT

- **Large selection of sizes and lengths** for the individual and perfect fit
- **Extensive range of accessories** for easy and convenient replacement of the gastrostomy probe – additional components not necessary
- **High comfort** due to soft, flexible silicone
- Low-profile design with optimised profile height, for **comprehensive mobility** and more quality of life
- **Easy cleaning** and maintenance

### SAFETY


- Made of **high-quality medical grade silicone** with high clarity
- All sets are **latex-free, PVC-free, DEHP-free, BPA-free**
- **Radiopaque stripe** along the entire shaft
- **Flat, rounded tip** enables atraumatic insertion
- **Recessed tip** to avoid trauma
- Secure anti-reflux valve
- Side valve access for balloon filling

The extension sets included in the danuButton® Compact and Premium Set are made of kink-resistant, biocompatible polyurethane. They are equipped with ENFit® connectors according to ISO 80369-3 and danumed's rotating connector, the ENSwivel®.



## The danuButton® Compact Set contains:


- danuButton® Gastrostomy Button
- danumed® Extension Set for Gastrostomy Button, 30 cm, Y-ENFit®, ENSwivel®, angled button connector
- danumed® Stiffener (Introducer)
- danumed® Balloon Syringe 10 ml
- optiLube lubricant
- Gauze Pads (x4)

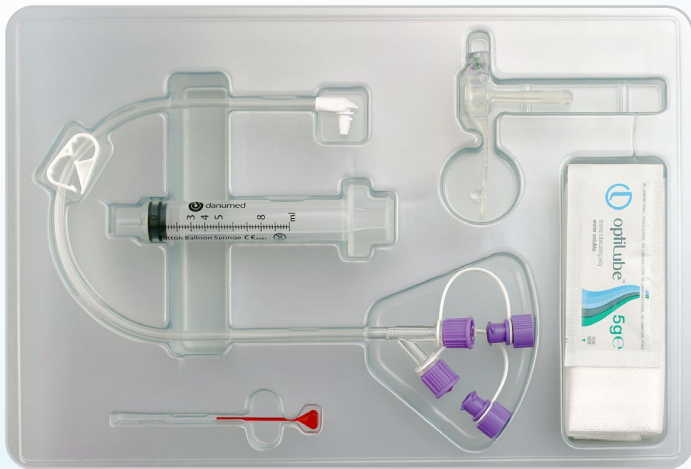

		Fr-Size					
		12	14	16	18	20	24
Length (mm)	8	DA 1042 1208	DA 1042 1408	DA 1042 1608	DA 1042 1808	DA 1042 2008	DA 1042 2408
	10	DA 1042 1210	DA 1042 1410	DA 1042 1610	DA 1042 1810	DA 1042 2010	DA 1042 2410
	12	DA 1042 1212	DA 1042 1412	DA 1042 1612	DA 1042 1812	DA 1042 2012	DA 1042 2412
	15	DA 1042 1215	DA 1042 1415	DA 1042 1615	DA 1042 1815	DA 1042 2015	DA 1042 2415
	17	DA 1042 1217	DA 1042 1417	DA 1042 1617	DA 1042 1817	DA 1042 2017	DA 1042 2417
	20	DA 1042 1220	DA 1042 1420	DA 1042 1620	DA 1042 1820	DA 1042 2020	DA 1042 2420
	23	DA 1042 1223	DA 1042 1423	DA 1042 1623	DA 1042 1823	DA 1042 2023	DA 1042 2423
	25	DA 1042 1225	DA 1042 1425	DA 1042 1625	DA 1042 1825	DA 1042 2025	DA 1042 2425
	27	DA 1042 1227	DA 1042 1427	DA 1042 1627	DA 1042 1827	DA 1042 2027	DA 1042 2427
	30	DA 1042 1230	DA 1042 1430	DA 1042 1630	DA 1042 1830	DA 1042 2030	DA 1042 2430
	35	DA 1042 1235	DA 1042 1435	DA 1042 1635	DA 1042 1835	DA 1042 2035	DA 1042 2435
40	DA 1042 1240	DA 1042 1440	DA 1042 1640	DA 1042 1840	DA 1042 2040	DA 1042 2440	
45	DA 1042 1245	DA 1042 1445	DA 1042 1645	DA 1042 1845	DA 1042 2045	DA 1042 2445	
50	DA 1042 1250	DA 1042 1450	DA 1042 1650	DA 1042 1850	DA 1042 2050	DA 1042 2450	

## The danuButton® Premium Set contains:

- danuButton® Gastrostomy Button
- danumed® Extension Set for Gastrostomy Button, 30 cm, Bolus-ENFit®, ENSwivel®, straight button connector
- danumed® Extension Set for Gastrostomy Button, 30 cm, Y-ENFit®, ENSwivel®, angled button connector
- danumed® Stiffener (Introducer)
- danumed® Balloon Syringe 10 ml
- danumed® Enteral Syringe 60 ml, ENFit®, female
- optiLube lubricant
- Gauze Pads (x4)

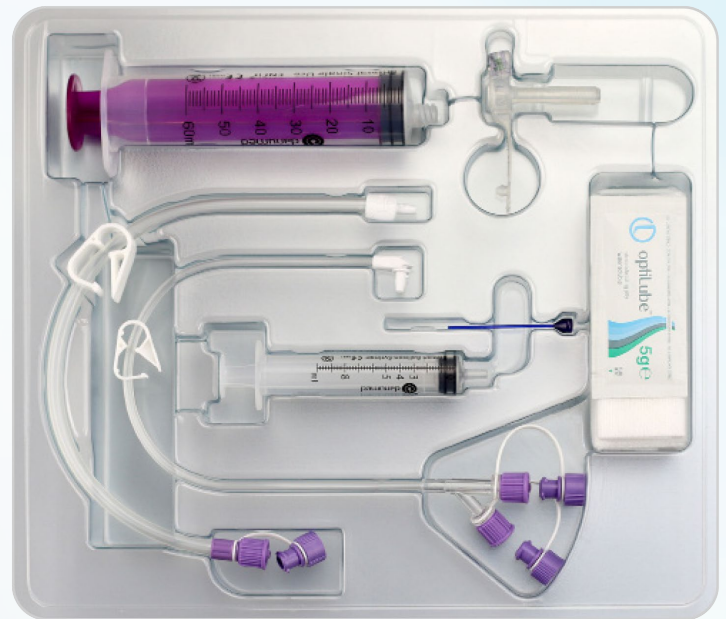

		Fr-Size					
		12	14	16	18	20	24
Length (mm)	8	DA 1041 1208	DA 1041 1408	DA 1041 1608	DA 1041 1808	DA 1041 2008	DA 1041 2408
	10	DA 1041 1210	DA 1041 1410	DA 1041 1610	DA 1041 1810	DA 1041 2010	DA 1041 2410
	12	DA 1041 1212	DA 1041 1412	DA 1041 1612	DA 1041 1812	DA 1041 2012	DA 1041 2412
	15	DA 1041 1215	DA 1041 1415	DA 1041 1615	DA 1041 1815	DA 1041 2015	DA 1041 2415
	17	DA 1041 1217	DA 1041 1417	DA 1041 1617	DA 1041 1817	DA 1041 2017	DA 1041 2417
	20	DA 1041 1220	DA 1041 1420	DA 1041 1620	DA 1041 1820	DA 1041 2020	DA 1041 2420
	23	DA 1041 1223	DA 1041 1423	DA 1041 1623	DA 1041 1823	DA 1041 2023	DA 1041 2423
	25	DA 1041 1225	DA 1041 1425	DA 1041 1625	DA 1041 1825	DA 1041 2025	DA 1041 2425
	27	DA 1041 1227	DA 1041 1427	DA 1041 1627	DA 1041 1827	DA 1041 2027	DA 1041 2427
	30	DA 1041 1230	DA 1041 1430	DA 1041 1630	DA 1041 1830	DA 1041 2030	DA 1041 2430
	35	DA 1041 1235	DA 1041 1435	DA 1041 1635	DA 1041 1835	DA 1041 2035	DA 1041 2435
40	DA 1041 1240	DA 1041 1440	DA 1041 1640	DA 1041 1840	DA 1041 2040	DA 1041 2440	
45	DA 1041 1245	DA 1041 1445	DA 1041 1645	DA 1041 1845	DA 1041 2045	DA 1041 2445	
50	DA 1041 1250	DA 1041 1450	DA 1041 1650	DA 1041 1850	DA 1041 2050	DA 1041 2450	

All danuButton® Sets are


- Latex-free • PVC-free • DEHP-free • BPA-free

# danumed® Extension Set for Gastrostomy Button

Extension set for delivery of enteral nutrition, liquid medication and decompression.

Innovative features optimise the flow rate and simplify the handling for the user.

Cost-efficient solution with usage for up to 14 days.

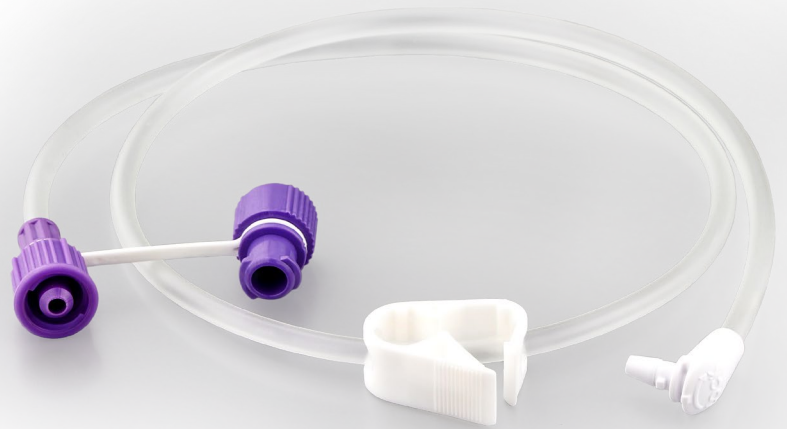

## COMFORT



### Rotating Connector



- Allows for the easy connection and release of ENFit® connections.
- Even very tight connections can be easily released
- Only the rotating Swivel is turned and not the entire connector. This feature **avoids any twisting of the enteral tube** which is both painful for the patient and inconvenient for the user
- All danumed ENFit® connectors offer the advantages of the ENSwivel®: **ENFit®**, **Y-ENFit®** and **Bolus-ENFit®**

### High-End Button Connector

- The locking mechanism with optimised design ensures **easy and tight connection** to the button

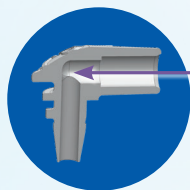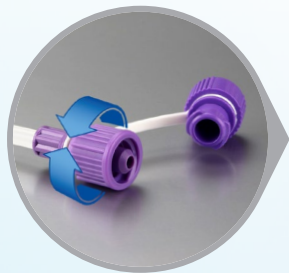

The rounded kink in the inner lumen optimises the flow rate, independent of the tube lumen.



Ergonomic form of the button connector with ribbed profile for a better grip.

## SAFETY

### Increased safety for the patient

- The ENFit® connector according to ISO 80369-3 **avoids tubing misconnections** and enhances patient safety significantly
- All ENFit® connectors from danumed are made of **TRITAN™**, a plastic material resistant to aggressive solutions/medicines
- All extension sets with ENFit® connector are **free from latex, PVC, DEHP and BPA**
- The biocompatible tube made of polyurethane is **kink resistant** and **reduces the hazard of tube occlusions** and the number of alarm signals of the enteral pump



## danumed® Extension Set with angled Button Connector



REF	Description	Length (cm)	Feeding connector	Latex-free	PVC-free	DEHP-free	BPA-free
DA 1011 3011	Extension Set, <b>sterile</b> , 30 cm, ENFit®, ENSwivel®	30 cm	ENFit®	●	●	●	●
DA 1011 6011	Extension Set <b>sterile</b> , 60 cm, ENFit®, ENSwivel®	60 cm	ENFit®	●	●	●	●
DA 1011 0521*	Extension Set <b>sterile</b> , 5 cm, Y-ENFit®, ENSwivel®	5 cm	Y-ENFit®	●	●	●	●
DA 1011 3021	Extension Set <b>sterile</b> , 30 cm, Y-ENFit®, ENSwivel®	30 cm	Y-ENFit®	●	●	●	●
DA 1011 6021	Extension Set <b>sterile</b> , 60 cm, Y-ENFit®, ENSwivel®	60 cm	Y-ENFit®	●	●	●	●

\*Item ordered on request

Delivered with right-angled button connector and **standard bore tube**. Individually packaged. 10 pieces in a box, 10 boxes in a shipping box.



## danumed® Extension Set with straight Button Connector



REF	Description	Length (cm)	Feeding connector	Latex-free	PVC-free	DEHP-free	BPA-frei
DA 1011 3032	Extension Set <b>sterile</b> , 30 cm, Bolus-ENFit®, ENSwivel®	30 cm	Bolus-ENFit®	●	●	●	●
DA 1011 6032	Extension Set <b>sterile</b> , 60 cm, Bolus-ENFit®, ENSwivel®	60 cm	Bolus-ENFit®	●	●	●	●

Delivered with straight button connector and **wide bore tube**. Individually packaged. 10 pieces in a box, 10 boxes in a shipping box. The Bolus Extension Sets are available on request with funnel connector instead of ENFit® connector.

### Related Products to danumed® Extension Sets:



danuButton®  
Gastrostomy Button

# danumed® GastroTube



## Gastrostomy probe with retention balloon for long-term enteral nutrition

*The danumed® GastroTube is a percutaneous gastric feeding and replacement tube made of silicone. It can replace a PEG tube after a completely healed stoma. Due to the wide range of sizes from 12 Fr to 24 Fr, almost every patient can be supplied with the perfectly fitting danumed® GastroTube.*

*Outstanding product features offer increased patient safety and best comfort in application and handling.*



## COMFORT

- **Extensive selection of sizes** from CH 12 to CH 24 for the individual and perfect fit
- **Patient comfort** due to soft, flexible silicone
- Soft and adjustable retention plate for **secure fixation**
- Structure of the retention plate allows better **air circulation** in the stoma site area
- **Length marking** at a distance of 2 cm
- **Display of the balloon volume** on the shaft
- **Side valve access** for balloon filling
- Valve colour coding for **easy sizing**
- **Convenient accessories** in every single package: a pack of **lubricant** for easy insertion and a **danumed® balloon syringe 10 ml** to block the retention balloon
- **Easy to clean** and maintain

## SAFETY

- Gastric replacement feeding tube made of **high-quality medical grade silicone**
- All danumed® GastroTubes are **free of Latex, PVC, DEHP and BPA**
- **Flat rounded tip** allows atraumatic insertion
- **Recessed tip** to prevent trauma
- High balloon stability of the inner retention balloon for **long usage and secure fixation**
- **ENFit®** connectors according to EN ISO 80369-3 avoid tubing misconnections with other systems
- Available with **ENFit®** and **Y-ENFit®** connectors
- **Radiopaque stripe** along the shaft

## danumed® GastroTubes with ENFit® and Y-ENFit®


REF	Product description	Size in Fr	Pieces/sales unit	Latex-free	PVC-free	DEHP-free	Bisphenol-A-free
DA 1044 0112	danumed® GastroTube with ENFit®	12	1 piece	●	●	●	●
DA 1044 0114	danumed® GastroTube with ENFit®	14	1 piece	●	●	●	●
DA 1044 0116	danumed® GastroTube with ENFit®	16	1 piece	●	●	●	●
DA 1044 0118	danumed® GastroTube with ENFit®	18	1 piece	●	●	●	●
DA 1044 0120	danumed® GastroTube with ENFit®	20	1 piece	●	●	●	●
DA 1044 0122	danumed® GastroTube with ENFit®	22	1 piece	●	●	●	●
DA 1044 0124	danumed® GastroTube with ENFit®	24	1 piece	●	●	●	●
DA 1044 0212	danumed® GastroTube with Y-ENFit®	12	1 piece	●	●	●	●
DA 1044 0214	danumed® GastroTube with Y-ENFit®	14	1 piece	●	●	●	●
DA 1044 0216	danumed® GastroTube with Y-ENFit®	16	1 piece	●	●	●	●
DA 1044 0218	danumed® GastroTube with Y-ENFit®	18	1 piece	●	●	●	●
DA 1044 0220	danumed® GastroTube with Y-ENFit®	20	1 piece	●	●	●	●
DA 1044 0222	danumed® GastroTube with Y-ENFit®	22	1 piece	●	●	●	●

Each danumed® GastroTube is individually sterile packed. The danumed® GastroTube is delivered in a sterile peel-pouch and additionally with lubricant and a sterile danumed® balloon syringe 10ml in each single box. 10 boxes in a shipping carton.

# danumed® StomaMeter

## Stoma Measuring Device with retention balloon

*The StomaMeter is a measuring instrument for determining the stoma length. The length of the stoma tract is measured simply, accurately and reliably to ensure the perfect selection of the optimal gastrostomy button.*

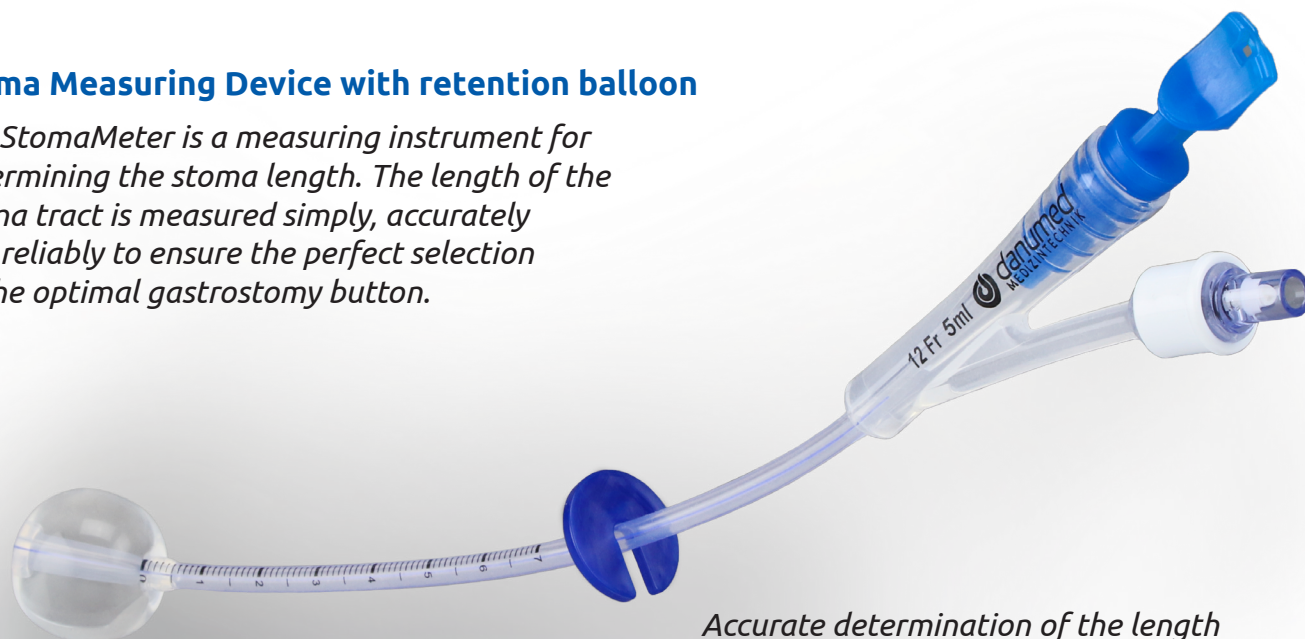

*Accurate determination of the length of the enteral stoma improves the fit of the gastrostomy button and thus increased patient safety and comfort.*

## COMFORT

- **Soft, flexible material** for best patient comfort
- **Non-slip retention plate** simplifies reading of measured values
- **Imprint of the balloon volume** on the funnel
- **Convenient accessories** in every single package: A pack of **lubricant** for easy insertion and a **danumed® balloon syringe 10 ml** to inflate the retention balloon.
- Available in **size 12 Fr**, suitable for each stoma



## danumed® StomaMeter

REF	Size in French	Pieces/sales unit
DA 1045 0112	12	1 piece

The danumed® StomaMeter stoma measuring device is individually sterile packed. The StomaMeter is delivered in a sterile peel-pouch and additionally with lubricant and a sterile danumed® balloon syringe 10ml in each single box. 10 boxes in a shipping carton.

The danumed® StomaMeter is

- **Latex-free**
- **PVC-free**
- **DEHP-free**
- **BPA-free**

## SAFETY

- Made of high-quality **medical grade silicone**
- **Latex-free, PVC-free, DEHP-free** and **BPA-free** material
- **Flat rounded tip** to avoid trauma
- **Radiopaque stripe** along the shaft
- **Precisely readable** and wipe-resistant scaling

For more information, please contact us on (02) 9882 3666 or at [amsl@amsl.com.au](mailto:amsl@amsl.com.au)

[amsl.com.au](http://amsl.com.au)

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE.  
AMSL is a Dexcom company. PR-100-741 July 2023

