


# Premium Plus™ Static Replacement Mattress

Product Code - APH001



The **Premium Plus** is a high specification replacement mattress designed to meet the demands of the acute healthcare environment. Practical and durable, this mattress delivers exceptional levels of comfort and pressure reduction. This is achieved by premium 'keyhole cut' foam that sits in a 'U' shaped core. The **Premium Plus** is suitable for those considered at **'High Risk'** of developing pressure ulcers.

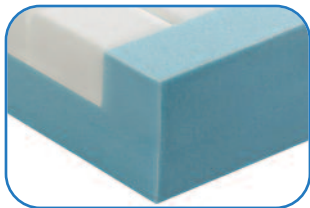

# Apollo Healthcare Technologies

- Featuring a unique contoured pad, constructed from 'keyhole cut' castellated foam cells, which move independently to provide excellent patient comfort, reduce shear and friction and optimise pressure reduction. In turn, this allows excellent conformity when profiled on electric bed frames.
- The outer cover is a multi-stretch, waterproof, vapour permeable, polyurethane fabric with welded seams and a full length flap concealing the zip. These features combine to effectively reduce the ingress of fluid into the mattress.
- The mattress cover can be wiped down or completely removed and laundered up to 95°C to aid decontamination protocols.

## Specifications

**Dimensions:** 198cm x 88cm x 15.6cm

**Maximum Weight Limit:** 39 stone / 247.6 Kg

**Risk Factor:** High Risk

**Warranty:** 5 years cover / 5 years foam

**Mattress Weight:** 14 Kg

**Safety Standards:** BS EN 597-1 Cig, BS EN 597-2 Match, BS7177: 1996 Crib Source 5



**SURGICAL HOUSE**

46 King Edward Road, Osborne Park WA 6017  
PO Box 1537, Osborne Park DC 6916

[surgical@surgicalhouse.com.au](mailto:surgical@surgicalhouse.com.au)  
[surgicalhouse.com.au](http://surgicalhouse.com.au)

+61 8 9382 3009  
+61 8 9381 4199