


# FINDING THE SUITABLE FIXOMULL® DRESSING

The comprehensive product range by **Fixomull®** offers **proven quality** for wide-area medical fixations. **Fixomull®** products are the right choice for the individual fixation needs of medical experts – **thanks to transparent, water-resistant and extremely gentle and skin-friendly dressings.** The entire range is available in rolls which can be cut to size and which are available in different widths and lengths for individual wound dressings.

## Fixomull® range of dressings

### Fixomull® stretch

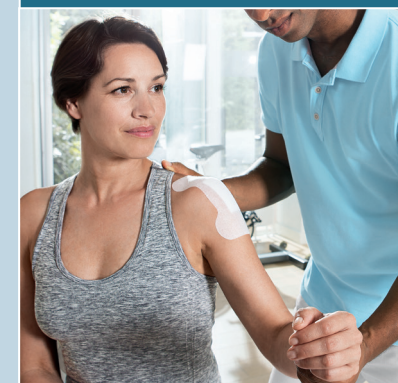

### Fixomull® transparent



### Fixomull® gentle touch



Indications and applications

- Transversely stretchable and particularly supple
- Especially suitable for joints and moving and contoured parts of the body
- For secure fixation of tubes, catheters, measuring devices and instruments

- Transparent, flexible and watertight
- Protects against moisture, germs and viruses
- Durable
- For secure fixation of tubes, catheters, measuring devices and instruments
- Hydrotherapy

- Particularly skin-friendly
- Atraumatic removal
- Easy repositioning
- Recommended for aging, cortisone, and parched skin and for neonatal care

- For particularly secure and wide-area fixation of wound dressings
- Good immediate and long-lasting adhesion
- Skin-friendly
- Easy application
- Economical in terms of materials and costs as it can be cut to size
- Radiotransparent

## ORDERING INFORMATION

Ref. No	Size	Order Qty	Shipping Unit
<b>Fixomull® gentle touch</b>			
79965-01	5cm x 5m	1 box	12 boxes
79965-02	10cm x 5m	1 box	12 boxes
79965-03	15cm x 5m	1 box	12 boxes
79965-04	10cm x 2m	1 box	12 boxes
<b>Fixomull® stretch</b>			
02033-00	15cm x 2m	1 box	12 boxes
02036-00	5cm x 10m	1 box	24 boxes
02037-00	10cm x 10m	1 box	24 boxes
02038-00	15cm x 10m	1 box	24 boxes
02039-00	20cm x 10m	1 box	24 boxes
02040-00	30cm x 10m	1 box	24 boxes
70022-01	10cm x 2m	1 box	12 boxes
<b>Fixomull® transparent</b>			
72216-00	5cm x 10m	1 box	12 boxes
72216-01	10cm x 10m	1 box	12 boxes
72217-00	10cm x 2m	1 box	12 boxes

Fixomull®

**FIXOMULL®**  
**GENTLE TOUCH**

SOFT AND SECURE FIXATION

BSN medical Australia & New Zealand

L6, 211 Wellington Road,  
Mulgrave, Victoria, Australia 3170

[www.bsnmedical.com.au](http://www.bsnmedical.com.au)  
[www.bsnmedical.co.nz](http://www.bsnmedical.co.nz)

BSN0543 01/2019

Customer Service Australia  
T 1300 276 633 F 1300 998 830

Customer Service New Zealand  
T 0508 276 111 F 0508 998 830

©Registered Trademark

**BSN**medical

an Essity company

**THERAPIES. HAND IN HAND.**

**BSN**medical

an Essity company

# FIXOMULL® GENTLE TOUCH – THE GENTLE AND SECURE FIXATION

We feel through our skin. Every touch, but every pain too. That is exactly why it is important to preserve and protect skin's integrity. When dealing with sensitive skin, greater attention is to be paid to withstand the stresses of everyday life.

To ease the care of patients with sensitive skin BSN medical has expanded the proven **Fixomull®** product range with the addition of a very gentle fixation bandage – **Fixomull® gentle touch**. Thanks to the **soft silicone adhesive** extremely sensitive skin such as baby's skin, aging skin, parchment skin or cortisone skin are well looked after. This way the skin is protected after and a successful treatment is provided.



## SKIN-FRIENDLY SILICONE ADHESIVE


Particularly supple on the skin surface

- Atraumatic removal
  - ▶ Repositioning without stressing the skin
  - ▶ Protects the skin and keeps the skin cells intact
- Secure fixation of dressings and splints
- Suitable as additional fixation of tubes, probes and catheters, amongst others
  - ▶ Supports the treatment process



## INDIVIDUAL SOLUTIONS FOR ALL WOUND SIZES

- **Fixomull® gentle touch** is available in rolls of varying widths
  - ▶ Can be cut to the desired size to treat any wound size
  - ▶ Quick and easy to apply
  - ▶ Economical and cost-effective

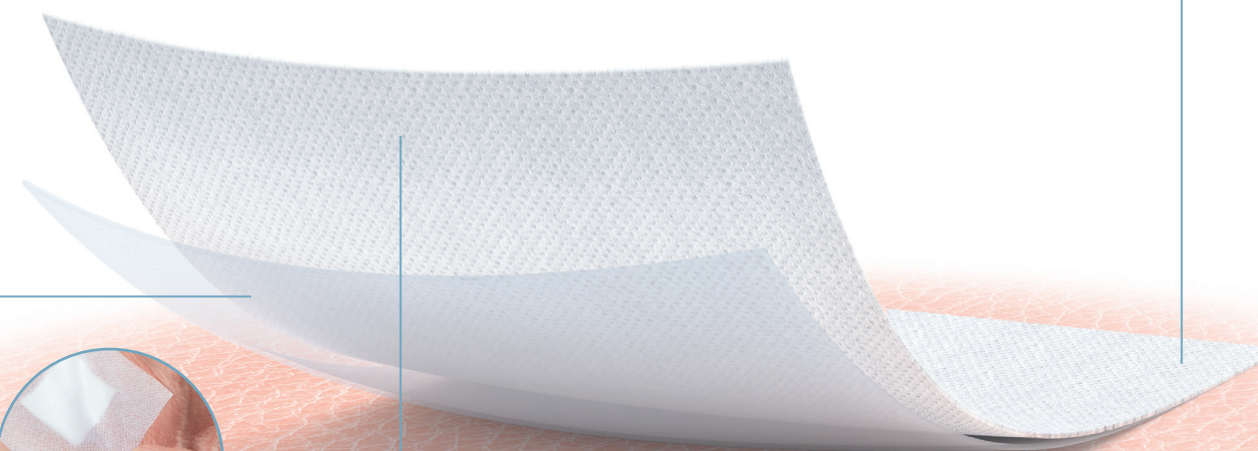

## SOFT NON-WOVEN MATERIAL

Highly flexible and transversely stretchable

- Can be applied to contoured parts of the body
  - ▶ Freedom of movement is barely restricted
  - ▶ Comfortable to wear
- Secure fixation of absorbent dressings

Breathable (water vapor permeable) surface

- ▶ Reduces the risk of skin macerations

Radiotransparent material

- Dressing does not have to be removed for X-ray examinations
  - ▶ Allows the wound to heal

# ALWAYS THE RIGHT CHOICE FOR SENSITIVE SKIN


## Indications

**Fixomull® gentle touch** is particularly suitable for sensitive skin:


- Aging skin
- Newborn/premature baby skin
- Parchment skin
- Dialysis patients
- Medicinally stressed skin

It can be applied to contoured parts of the body and can be repositioned several times. **Fixomull® gentle touch** is recommended for the secure fixation of dressings and the secondary fixation of probes, tubes, measuring devices and instruments. It can be worn for up to seven days.


## APPLICATION INSTRUCTIONS FOR FIXOMULL® GENTLE TOUCH

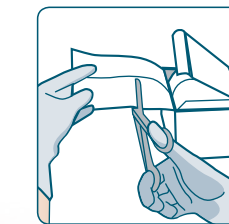

Cut a piece of **Fixomull® gentle touch** to the required size.

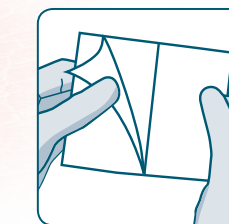

The backing sheet is separated. Remove the first part.

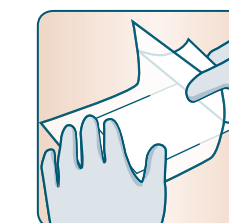

Place the bandage over the dressing and carefully apply to the skin. Fold the other side back and place the rest of the bandage over the dressing.

### Wound specialist tip

**Fixomull® gentle touch** can be used as a fenestrated dressing in order to securely fix non-adhesive dressings.

