



ASPIRE STAR PATIENT MOVER

The Aspire Star Mover offers a safe and comfortable solution for patients and carers doing sit to stand transfers. The compact and versatile design is perfect for Homecare, Institutional and Aged Care environments.

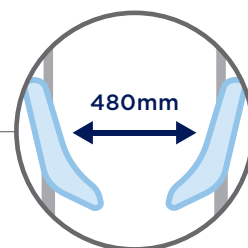

Wide Seat Entry Space

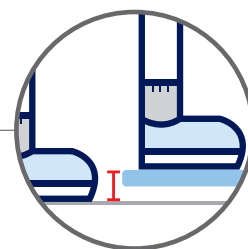

Low Footplate Height



Large Castors



2 YEAR WARRANTY



OVERALL LENGTH



OVERALL HEIGHT



OVERALL WIDTH



KNEE PAD HEIGHT



SEAT PAD HEIGHT



HEIGHT HANDLE



LEG RANGE



FOOTPLATE HEIGHT FROM FLOOR

LSS390580 Star Mover	925mm	1030mm	615mm - 900mm	205mm - 495mm	685mm - 850mm	1030mm	Min External Leg Width: 520mm Max Internal Leg Width: 740mm	60mm
	<small>*Range of motion dimensions are accurate to a +/-10mm tolerance</small>							

FEATURES & BENEFITS

- Wide seat entry suits a large range of body shapes
- Low footplate height makes patient entry and exit easier
- Large 100mm castors (two rear locking) have low rolling resistance and offer a smooth ride during transportation
- Narrow base improves access in tight spaces
- Height adjustable knee support ensures correct positioning for a larger patient height range
- Easy access seat pad handles provide reach-around access to carers
- Soft touch grab bar provides a comfortable patient support for improved stability and confidence during transfers
- Width adjustable legs with tapered close position allow access between and, or around chair legs

HOW TO USE



Open legs if required and lift seat pads. Ask the user to place feet on footplate and apply brakes.




Ask the user to hold the grab handle and stand up.



Swing the seat pads back into place and ask the user to sit down. Unlock brakes and commence the transfer.

SLING COMPATIBILITY

The Star Mover features clip attachment points that work with Aspire Deluxe Standing Slings.

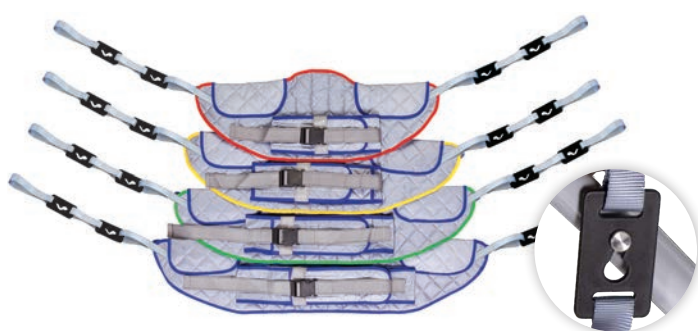


Code	Description	Girth	Edging	S.W.L
LSS390920	Deluxe Padded Standing Sling - Polyester - Clip - Small	65-95cm	Red	200kg
LSS390925	Deluxe Padded Standing Sling - Polyester - Clip - Medium	77-110cm	Yellow	200kg
LSS390930	Deluxe Padded Standing Sling - Polyester - Clip - Large	98-135cm	Green	300kg
LSS390935	Deluxe Padded Standing Sling - Polyester - Clip - Extra Large	118-157cm	Blue	300kg

For men and children, girth is the chest circumference. For women, girth is the hip circumference.



Star Patient Mover is registered in the ARTG, AUST R ARTG 373695

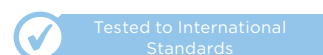

Tested to International Standards

ASFY71.6 221021