

General information and sizing

All LiftAssist® slings are tested by an accredited NATA Laboratory to meet Australian Standards AS/NZ ISO 10535:2011

with a safe working load of:

Sizes:


XXS-L 250 kgs / 550 lbs
 XL - XXL 363 kgs/ 800 lbs

TGA registration is ARTG 235001

The likely size of sling to meet the client's needs and the colour coded tape/webbing which indicates sling size are shown on the chart (see right).

Size	Size Code	Height	Weight
XXS	01	< 116 cm < 3'10"	< 40 kg < 88 lbs
Junior	03	< 121 cm < 4'	< 49 kg < 110 lbs
Small	06	121 cm - 167 cm 4' - 5'6"	43 kg - 68 kg 95 lbs - 150 lbs
Medium	09	152 cm - 182 cm 5' - 6'	56 kg - 113 kg 125 lbs - 250 lbs
Large	12	182 cm - 213 cm 6' - 7'	113 kg - 181 kg 250 lbs - 400 lbs
X Large	15	182 cm - 213 cm 6' - 7'	181 kg - 217 kg 400 lbs - 480 lbs
XX Large	18	> 198 cm >6'6"	> 217 kg >480 lbs

Comfort Recline

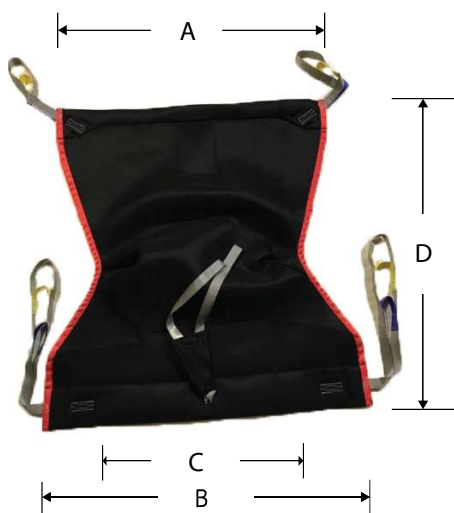

The LiftAssist® Comfort Recline is an 'ALL DAY' sling made possible by the absence of seams and the use of Spacer fabric which has superior stretch and breathing properties enabling it to be left in place under the person. Other types of fabric increase the risk of skin breakdown if left underneath the person for long periods of time. The Comfort Recline sling can only be applied while the person is supine, it is mainly used for supine to seated transfers, remains in place for long periods, then transfer from seated back to supine.

Size	Colour Code	A	B	C	D	Item code
Junior	Grey	690	740	510	620	701203
Small	Red	770	810	590	850	701206
Medium	Yellow	890	905	665	940	701209
Large	Green	1050	1095	825	1075	701212
X Large	Blue	1180	1260	1010	1310	701215

Comfort Recline with head support



The LiftAssist® Comfort Recline with head support is an 'ALL DAY' sling made possible by the absence of seams and the use of Spacer fabric which has superior stretch and breathing properties enabling it to be left in place under the person. Other types of fabric increase the risk of skin breakdown if left underneath the person for long periods of time. The Comfort Recline sling can only be applied while the person is supine, it is mainly used for supine to seated transfers, remains in place for long periods, then transfer from seated back to supine.

Size	Colour Code	A	B	C	D	Item code
XXS	Orange	640	720	480	760	701301
Junior	Grey	685	740	510	835	701303
Small	Red	760	810	590	1140	701306
Medium	Yellow	855	905	665	1235	701309
Large	Green	825	1065	825	1390	701312
X Large	Blue	1180	1260	1010	1520	701315

Hygiene Sling



The LiftAssist® Hygiene features a U-shaped design to allow access to remove and replace clothing for toileting while still allowing for support for the back. The padded Tri-cot fabric makes for a more comfortable transfer for the client.

Size	Colour Code	A	B	C	D	Item code
Junior	Grey	840	265	495	760	700503
Small	Red	910	290	545	835	700506
Medium	Yellow	985	320	595	915	700509
Large	Green	1065	345	645	985	700512
X Large	Blue	1140	370	700	1070	700515
XX Large	Black	1220	395	450	1145	700518

General Purpose Hygiene with Head support



The LiftAssist® Gen Purpose Hygiene sling features a U-shaped design to allow access to remove and replace clothing for toileting. Our Hygiene sling with Head support adds that extra support for those with poor head or neck control. The body is made with mesh fabric which allows for air flow and quick drying times, but the padded legs which provide extra comfort requires extra drying time.

Size	Colour Code	A	B	C	D	Item code
Junior	Grey	860	540	495	1035	700703
Small	Red	940	610	545	1155	700706
Medium	Yellow	1010	690	595	1285	700709
Large	Green	1100	790	645	1435	700712
X Large	Blue	1200	880	700	1580	700715
XX Large	Black	1300	980	750	1730	700718


Universal Padded



The LiftAssist® Universal padded uses Tricot fabric and with stretch padding support makes it the most comfortable sling for general transfer. It is suitable for a wide variety of applications. Its split leg enables the users legs to be both more open or closed with the slings under each leg or positioned under both as a hammock.

Size	Colour Code	A	B	C	D	Item code
Junior	Grey	760	760	470	995	700303
Small	Red	830	820	550	1135	700306
Medium	Yellow	950	920	630	1275	700309
Large	Green	1100	1020	710	1420	700312
X Large	Blue	1280	1120	800	1575	700315
XX Large	Black	1400	1185	890	1715	700318


Universal padded with head support

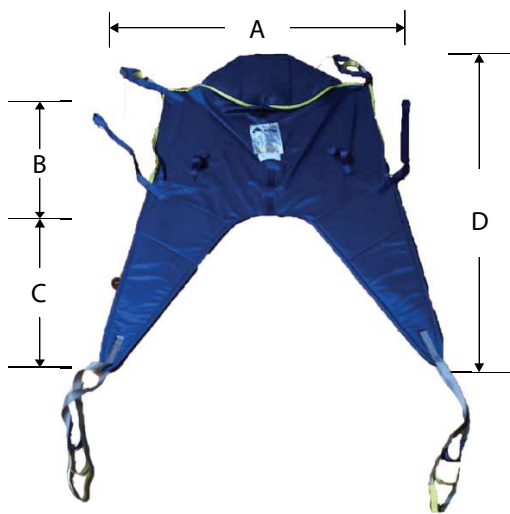

The LiftAssist® Universal padded with head support for much more on head control uses a Tricot fabric and with stretch padding support makes it the most comfortable sling for general transfer. It is suitable for a wide variety of applications. Its split leg enables the users legs to be both more open or closed with the slings under each leg or positioned under both as a hammock.

Size	Colour Code	A	B	C	D	Item code
Junior	Grey	760	525	470	1320	700403
Small	Red	890	585	550	1410	700406
Medium	Yellow	1030	645	630	1560	700409
Large	Green	1160	710	710	1740	700412
X Large	Blue	1280	775	800	1960	700415
XX Large	Black	1400	825	890	2195	700418

Universal mesh

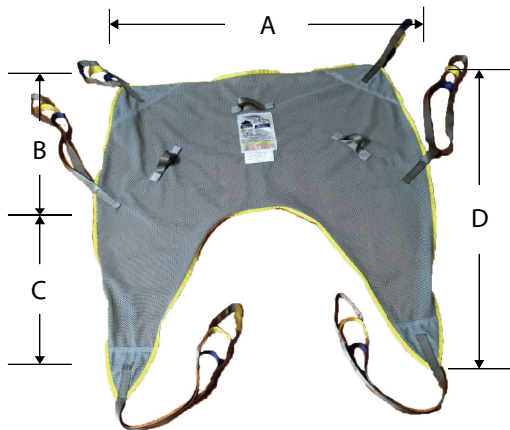

The LiftAssist® Universal sling is suitable for a wide variety of applications. Its split leg enables the user's legs to be both more open or closed with the slings under each leg or positioned under both as a hammock. It is most commonly supplied for general transfer purposes. Fully in mesh it is perfect for quick drying.

Size	Colour Code	A	B	C	D	Item code
Junior	Grey	760	525	430	955	700103
Small	Red	890	585	510	1095	700106
Medium	Yellow	1030	645	585	1230	700109
Large	Green	1160	710	710	1420	700112
X Large	Blue	1280	775	840	1615	700115
XX Large	Black	1400	825	965	1790	700118

Universal mesh with head support



The LiftAssist® Universal sling with head support for weak neck or head control is suitable for a wide variety of applications. Its split leg enables the user's legs to be both more open or closed with the slings under each leg or positioned under both as a hammock. It is most commonly supplied for general transfer purposes. Fully in mesh it is perfect for quick drying.

Size	Colour Code	A	B	C	D	Item code
Junior	Grey	760	525	430	1185	700203
Small	Red	890	585	510	1325	700206
Medium	Yellow	1030	645	585	1460	700209
Large	Green	1160	710	710	1650	700212
X Large	Blue	1280	775	840	1845	700215
XX Large	Black	1400	825	965	2020	700218


Universal General purpose deluxe mesh with head support



The LiftAssist® Universal GP deluxe Mesh with head support is suitable for a wide variety of applications. Its split leg enables the user's legs to be both more open or closed with the slings under each leg or positioned under both as a hammock. The body is mesh which enables air flow and quick drying but the padding in the legs and head support provides greater comfort and support but needs more drying time.

Size	Colour Code	A	B	C	D	Item code
Junior	Grey	760	760	430	1190	700803
Small	Red	830	820	510	1330	700806
Medium	Yellow	950	920	585	1505	700809
Large	Green	1100	1020	710	1730	700812
X Large	Blue	1280	1120	840	1960	700815
XX Large	Black	1400	1185	910	2095	700818

Hammock Mesh

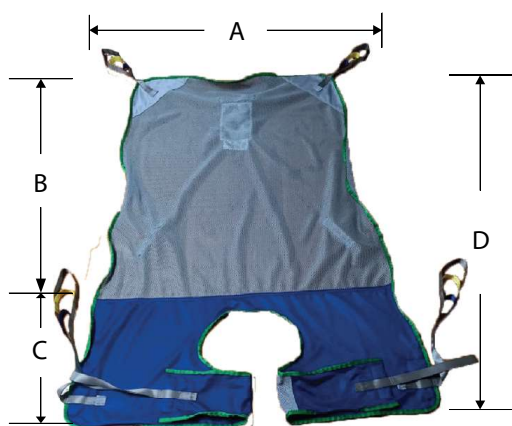

The LiftAssist® Hammock Mesh is designed for patients that have poor or weakened body tone. It provides maximum support to the sacral area and is used in configurations of divided leg open, divided leg crossed or closed leg. This sling also offers appropriate support for amputees.

Size	Colour Code	A	B	C	D	Item code
Junior	Grey	635	455	355	815	700603
Small	Red	685	635	380	1015	700606
Medium	Yellow	735	815	405	1220	700609
Large	Green	785	990	405	1400	700612
X Large	Blue	840	1170	405	1575	700615
XX Large	Black	890	1345	405	1755	700618


Hammock Spacer



This LiftAssist® Hammock Sling is made from Spacer fabric which has superior stretch and breathing properties enabling it to be left in place under the person. Other types of fabric increase the risk of skin breakdown if left underneath the person for long periods of time. It has a commode hole for toileting and split leg for different support methods.

Size	Colour Code	A	B	Item code
Junior	Grey	700	1040	701403
Small	Red	810	1090	701406
Medium	Yellow	930	1150	701409
Large	Green	1110	1220	701412
X Large	Blue	1210	1290	701415


Walking Sling

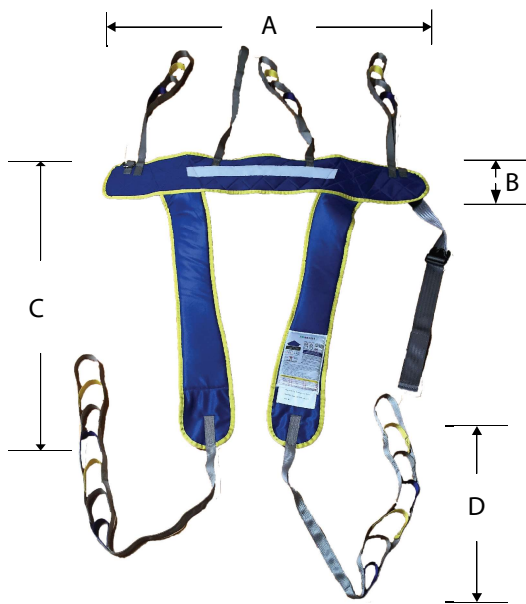

The LiftAssist® Walking Sling is designed to allow a person to be lifted vertically from a sitting position and positioned upright for walking. It is well suited for use with an assistive walking device such as parallel bars, hand rails etc.

Size	Colour Code	A	B	C	D	Item code
Medium	Blue	1016	177	1041	939	701110


Positioning Sling

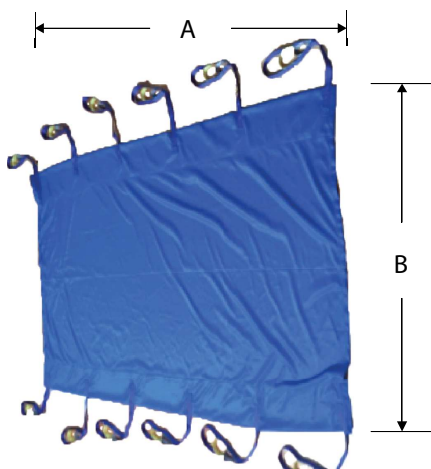

The LiftAssist® Positioning sling is designed for ceiling hoists with 2, 4, or 6 point standard carry bars with webbing loop attachment. It can remain on the bed for continuous use and can be used for repositioning in bed, transfer to a flat surface and for turning.

Size	Colour Code	A	B	Item code
Standard	Blue	1950	1050	701010


Tri Turner Sling

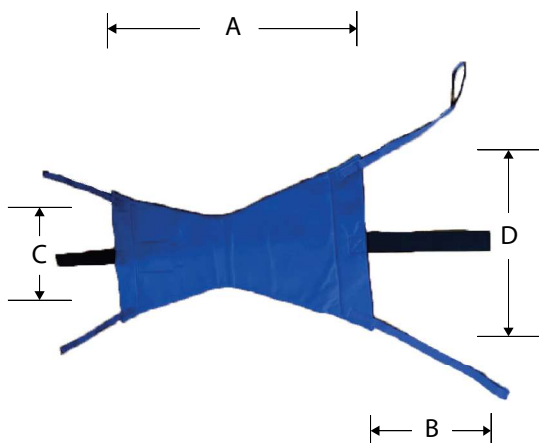

The LiftAssist® Tri Turner sling is to be used for turning purposes only and not for lifting. It has a unique design with a Velcro loop and hook to attach to bed frame on the opposite side to the hoist. The webbing loops attach to your carry bar.

Size	Colour Code	A	B	C	D	Item code
Medium	Blue	990	609	177	482	700909

Sit to Stand Sling

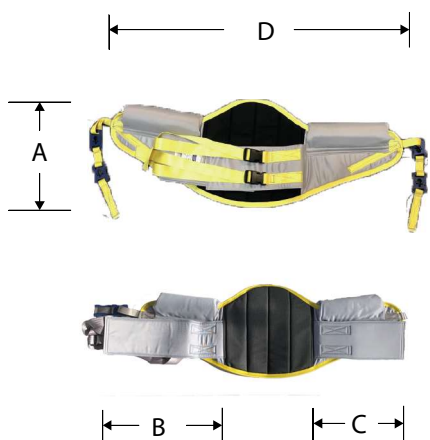

This LiftAssist® Sling is designed for Sit to Stand hoists that use compatible size “keyhole clips”. Manual transportation devices will require greater ability for the person being moved to stand with minimal or no assistance, even though they may be at risk of falling if walking.


This is only suitable for people who can safely weight-bear, but may be unsteady walking.

Size	Colour Code	A	B	C	D	Item code
Small	Red	310	370	270	880	701506
Medium	Yellow	385	490	370	1020	701509
Large	Green	395	620	520	1260	701512

Universal Sling with Head support and clips



The LiftAssist® Universal sling with clips is designed for hoists using a pivot / wishbone / DPS style carry bar with compatible “key hole” clip attachments. It has head support for those who have weak or limited neck and head control. It is available in 3 fabrics: Mesh - for quick drying; Tricot - for general comfort; Spacer - for its stretch comfort and breathability to assist leaving in place for lengthy periods.

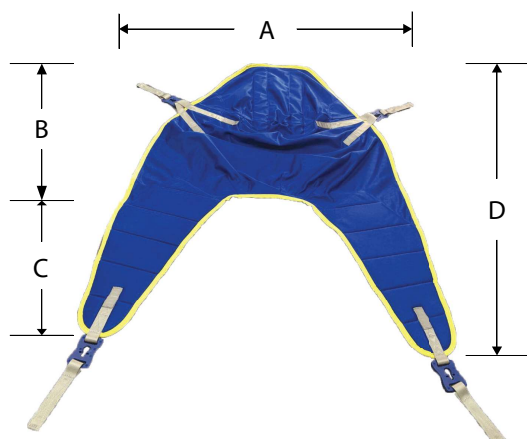

Size	Colour Code	A	B	C	D	Item code		
						Mesh	Tricot	Spacer
Small	Red	780	860	450	1130	701606	701706	701806
Medium	Yellow	790	890	540	1260	701609	701709	701809
Large	Green	940	930	600	1370	701612	701712	701812
X Large	Blue	1020	1050	625	1440	701615	701715	701815

Rehab Assist
2/3 Kearney Street
Bayswater, VIC, 3153
Ph: (03) 9729 7978
Email: sales@rehabassist.com.au
Website: www.rehabassist.com.au

Rehab Assist
Making Life Better